

LITHOSTAR MODULARIS

SP

Replacement instructions

Labeling the Frame Cover on the SONOCUR and
LITHOSTAR MODULARIS

© Siemens AG 2000

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

English

Doc. Gen. Date: 05.00

Chapter	Page	Revision
0	all	01
1	all	01
2	all	01
3	all	01

	Page
1 _____ General information _____	1 - 1
Tools.	1 - 1
2 _____ SONOCUR _____	2 - 1
Replacing the covers	2 - 1
3 _____ LITHOSTAR MODULARIS _____	3 - 1
Replacing the covers	3 - 1

This page intentionally left blank.

NOTICE	This kit contains the necessary labels for the SONOCUR and LITHOSTAR MODULARIS. Any unused labels should be destroyed.
---------------	---

- Contents of label kit 47 83 226

Label	Quan- tity	SONOCUR	LITHOSTAR MODULARIS
"Voltage indica- tor" sign	1	beside the power cable inlet	beside the power cable inlet
"Hospital grade" sign	1	-----	beside the power cable inlet
Revision status	1	beside the power cable inlet	beside the power cable inlet

Tools

- Protective conductor measuring device

This page intentionally left blank.

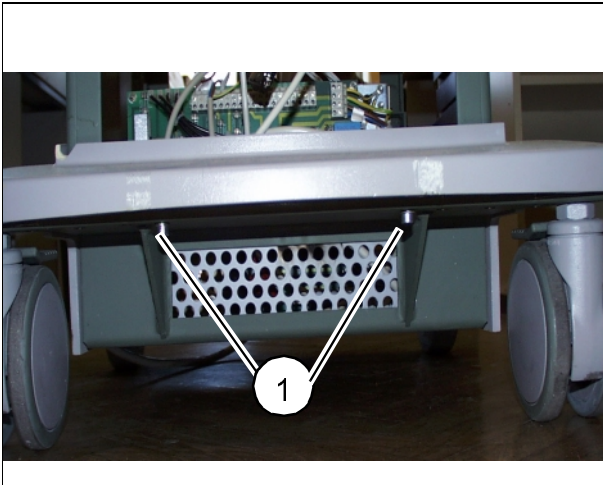


Fig. 1

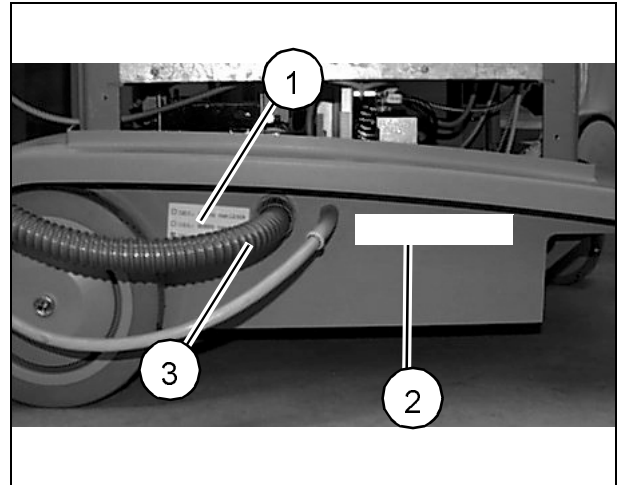


Fig. 2

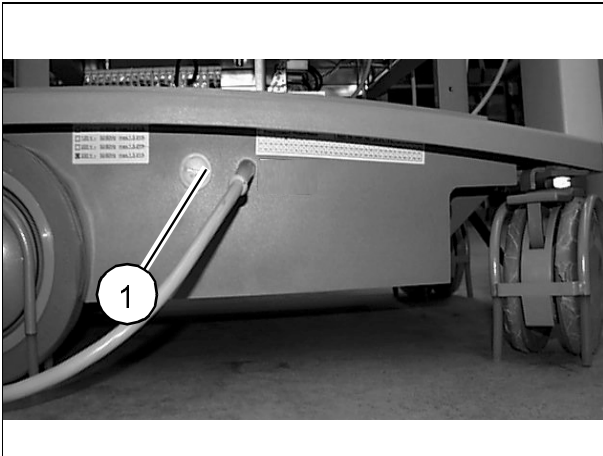


Fig. 3

Replacing the covers

- Remove the white cover completely.
- Unscrew the blue grip on the trolley.
- Remove the power cable.
- If the ultrasound option is available:
 - remove the grey hose (3/ Fig. 2)
- Loosen the covers (1/ Fig. 1); two screws are also located on the other side.
- Remove the covers by lifting them up (no component of the suspension arm needs to be removed).
- Reverse the order of these steps to attach the new covers.
- Check off the corresponding line voltage on the line voltage label (1/ Fig. 2) by referring to the old label and affix to the covers.

- Copy the information on the old label (model number and revision status) to the new revision status label (2/Fig. 2) and affix to the new covers.
- If the ultrasound option is not available:
 - Remove the steel cover (1/Fig. 3) from the old covers and attach it to the new covers.
- Check the protective conductor.

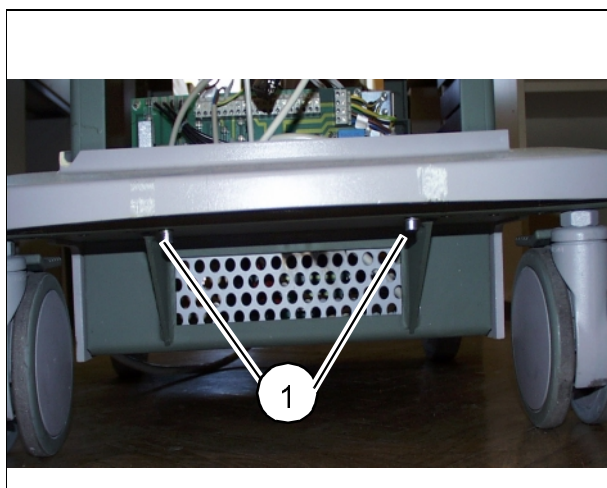


Fig. 1

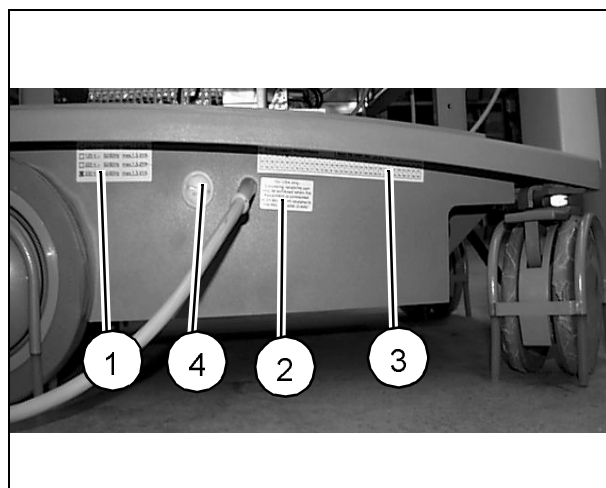


Fig. 2

Replacing the covers

- Remove the white cover completely.
- Unscrew the blue grip on the trolley.
- Remove the power cable.
- Loosen the covers (1/Fig. 1); two screws are also located on the other side.
- Remove the covers by lifting them up (no component of the suspension arm needs to be removed).
- Reverse the order of these steps to attach the new covers.
- Check off the corresponding line voltage on the line voltage label (1/Fig. 2) by referring to the old label and affix to the covers.
- Copy the information on the old label (model number and revision status) to the new revision status label (2/Fig. 2) and affix to the new covers.
- Affix the "Hospital grade" label (2/Fig. 2) to the covers.
- Remove the steel cover (4/Fig. 2) from the old covers and attach it to the new covers.
- Check the protective conductor.

This page intentionally left blank.